BARD® VENTRALEX™ Hernia Patch

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DAVOL INC.

BARD® VENTRALEX™ Hernia Patch

Self-Expanding Polypropylene & ePTFE Patch
with Strap and SORBAFLEX™ Memory Technology
for Soft Tissue Reconstruction

DESCRIPTION:


The BARD® VENTRALEX™ Hernia Patch is a self-expanding, nonabsorbable, sterile prosthesis. Two layers of monofilament polypropylene mesh are stitched with PTFE monofilament to an ePTFE sheet. The top layer of the polypropylene mesh is designed to form a strap and pocket. The strap facilitates placement, positioning and fixation of the device. The pocket is designed to aid in the fixation of the patch to the abdominal wall. Polypropylene monofilaments are knitted into a mesh, which allows a prompt fibroblastic response through the interstices of the mesh. The ePTFE layer has a sub-micronal porosity, which helps minimize tissue attachment. A depth marker on the positioning strap for the small BARD® VENTRALEX™ Hernia Patch (4.3 cm / 1.7") is designed to facilitate placement of the small patch through a trocar. The device contains SORBAFLEX™ Memory Technology which provides memory and stability to the device, facilitating ease of initial insertion, proper positioning, and fixation of the device. The SORBAFLEX™ Memory Technology is comprised of an extruded polydioxanone (PDO) absorbable monofilament that is contained within a knitted polypropylene mesh tube. The SORBAFLEX™ PDO monofilament fully degrades in vivo by means of hydrolysis. The PDO monofilament has been found to elicit an inflammatory response during absorption. Absorption is essentially complete in 6-8 months. The PDO monofilament is dyed violet by adding D & C Violet No. 2.

INDICATIONS:

The BARD® VENTRALEX™ Hernia Patch is intended for use in all forms of hernia repair requiring reinforcement with a nonabsorbable support material. The small BARD® VENTRALEX™ Hernia Patch (4.3 cm / 1.7") is also intended to repair soft tissue deficiencies, including deficiencies caused by trocars.

CONTRAINDICATIONS:

1. Do not use the BARD® VENTRALEX™ Hernia Patch in infants or children, whereby future growth will be compromised by use of such mesh material.
2. Do not use the BARD® VENTRALEX™ Hernia Patch for the reconstruction of cardiovascular defects.
3. Literature reports that there is a possibility for adhesion formation when the polypropylene is placed in contact with the bowel or viscera.

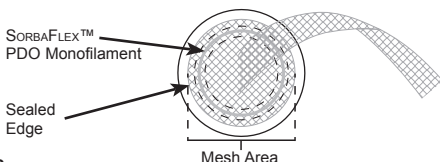WARNINGS:

1. The device is supplied sterile. Inspect the packaging to be sure it is intact and undamaged prior to use.
2. This device is for single use only. Do not resterilize or reuse any portion of the BARD® VENTRALEX™ Hernia Patch.
3. Do not cut or reshape any portion of the BARD® VENTRALEX™ Hernia Patch (as this could affect its effectiveness), except for the monofilament polypropylene positioning strap. Care should be taken not to cut or nick the SORBAFLEX™ PDO Monofilament. If the recoil ring is cut or damaged during insertion or fixation, additional complications may include bowel or skin perforation and infection.
4. Follow proper rolling techniques for all patches as described in these instructions for use as other rolling techniques may potentially compromise the SORBAFLEX™ PDO Monofilament.
5. Ensure proper orientation; the solid white surface (ePTFE) must be oriented against the bowel or sensitive organs. Do not place the mesh surface against the bowel. There is a possibility for adhesion formation when mesh (including strap) is placed in direct contact with the bowel or viscera.
6. To ensure a strong repair, the prosthesis should be secured with tacks or sutures through the polypropylene mesh straps and/or positioning pocket. Suturing or tacking on the sealed edge of mesh alone is not recommended.
7. Excess positioning strap material above the fixation line and at the level of the fascia must be cut off and discarded to eliminate excess material from remaining in the body.
8. When used to repair deficiencies caused by trocars, the device should be used under endoscopic guidance or direct visualization.
9. The use of any permanent mesh or patch in a contaminated or infected wound could lead to fistula formation and/or extrusion of the prosthesis.
10. If an infection develops, treat the infection aggressively. Consideration should be given regarding the need to remove the patch. An unresolved infection may require removal of the device.
11. To prevent recurrences when repairing hernias, the prosthesis should be large enough to extend beyond the margins of the defect.

PRECAUTIONS:

1. Please read all instructions prior to use.

- Only physicians qualified in the appropriate surgical techniques and use of this device should use this prosthesis.
- Care should be taken not to cut or nick the SORBAFLEX™ PDO Monofilament or the knitted polypropylene mesh tube during fixation.
- The safety and effectiveness of BARD® VENTRALEX™ Hernia Patch has not been evaluated in clinical studies in the presence of malignancies in the abdominopelvic cavity.

ADVERSE REACTIONS:

Possible complications include seroma, adhesions, hematoma, inflammation, extrusion, fistula formation, infection, allergic reaction, and recurrence of the hernia or soft tissue defect.

If the SORBAFLEX™ PDO Monofilament is cut or damaged during insertion or fixation, additional complications may include bowel or skin perforation and infection.


INSTRUCTIONS FOR USE:

Surface Orientation:

It is extremely important that this product is oriented correctly to function as intended. The smooth, solid white (ePTFE) surface is designed with low porosity to minimize tissue attachment. Place this side of the prosthesis against those surfaces where minimal tissue attachment is desired, i.e., against bowel or other visceral structures. It is still recommended to pull down omentum wherever possible beneath the ePTFE side to further mitigate the risk of visceral adhesion. The porous mesh side offers the same tissue ingrowth characteristics of polypropylene mesh alone. Therefore, this surface should face the surface where tissue ingrowth is desired. The mesh surface should never be placed against the bowel or other visceral structures.


Hernia Repair:

Select the patch that is approximately twice the size of the defect. Roll the patch parallel to the opening between the straps with the ePTFE side facing out for insertion into the defect. The free or unattached end of the positioning strap is left outside of the body during placement. The positioning strap is manipulated to facilitate the proper positioning of the patch. Gently pull up on the positioning strap, which will allow the patch to flatten itself against the underside of the abdominal wall. The mesh positioning strap can be pulled apart to gain access to the inner positioning pocket.

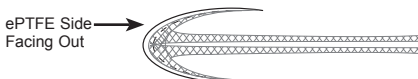Fixation:

To secure the patch, the mesh positioning strap is pulled up and apart and the device is then secured to the margins of the defect through the strap or positioning pocket. Bard permanent or absorbable fixation devices or nonabsorbable monofilament sutures are recommended to properly secure the prosthesis. If other fixation devices are used, they must be indicated for use in hernia repair. Care should be taken to ensure that the patch is adequately fixated to the abdominal wall. If necessary additional fasteners and/or sutures should be used. Excess positioning strap material above the fixation line and at the level of the fascia must be cut off and discarded. The method of securing the prosthesis should be determined by surgeon preference and the nature of the reconstruction to provide for adequate tissue fixation and to prevent reherniation.

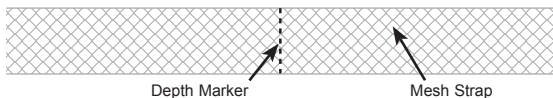Trocar Defect Repair:

Roll the small BARD® VENTRALEX™ Hernia Patch (4.3 cm / 1.7") with the ePTFE side facing out and hold it in the rolled position with an endoscopic grasper. Deploy the prosthesis all the way through the trocar and into the intraabdominal space. For a reference point, the product is constructed with a depth marker (blue line) on the positioning strap 8 inches (203 mm) above the patch. The size of the patch may inhibit deployment through small trocars. A minimum 10 mm diameter trocar should be used to introduce the patch. Release the prosthesis from the grasper when the patch is completely through the trocar. The SORBAFLEX™ PDO Monofilament allows the patch to open and lay flat in the intra-abdominal space. Pull up on the positioning strap to hoist the prosthesis up to the distal end of the trocar. Pull out the trocar and pass it over the strap. Pull up firmly without creating tension, on the positioning strap to make sure that the prosthesis is completely covering the defect. While pulling up on the positioning strap, suture close the anterior fascia, catching the positioning strap in between the margins of the fascia. Cut off excess positioning strap immediately above the suture at the fascia level and discard.

TRACEABILITY:

A traceability label which identifies the type, size and lot number of the prosthesis is attached to every package. This label should be affixed to the patient's permanent medical record to clearly identify the device which was implanted.

Bard, Davol, and Ventralex are trademarks and/or registered trademarks of C. R. Bard, Inc. or an affiliate.

Copyright © 2005-2008, 2011, 2013 C. R. Bard, Inc. All Rights Reserved.

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Copyright © 2005-2008, 2011, 2013 C. R. Bard, Inc. Tous droits réservés.
Copyright © 2005–2008, 2011, 2013 C. R. Bard, Inc. Alle Rechte vorbehalten.
Copyright © 2005-2008, 2011, 2013 C. R. Bard, Inc. Tutti i diritti riservati.
Copyright © 2005-2008, 2011, 2013 C. R. Bard, Inc. Todos los derechos reservados.
Copyright © 2005-2008, 2011, 2013 C. R. Bard, Inc. Alle rechten voorbehouden.
Copyright © 2005-2008, 2011, 2013 C. R. Bard, Inc. Todos os direitos reservados.