
Puritan Bennett™

**Re/X700 Reusable Expiratory
Filter for 700 Series Ventilators**

en	Instructions For Use
fr	Mode d'emploi
de	Gebrauchsanweisung
nl	Gebruiksaanwijzing
it	Istruzioni per l'uso
es	Instrucciones de uso
pt	Instruções de uso
ru	Инструкция по применению
pl	Instrukcja użytkowania

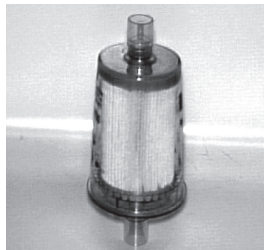Re/X700

Reusable Expiratory Filter for 700 Series Ventilators

REF # G-060525-00

 Indicates Flow direction. The Re/X700 is designed to be used only in this direction.
 Anticipated filter retirement date to be inserted here.

Directions for use

Description

The Re/X700 reusable expiratory filter reduces contaminants delivered to the patient, and minimizes the possibility of contaminants re-entering the ventilator, where they can cause reinfection or cross-infection. The filter is designed for repeated autoclaving. The see-through plastic shell allows visual inspection of the filter.

Specifications

Filter medium: Minimum filter efficiency of 99.97% retention of particles 0.3 μm nominal at 100 L/min flow.

Resistance to air flow: Less than 2.5 cmH_2O at 100 L/min flow when new.

Size: 3.5 in. diameter (8.9 cm) x 6.70 in. length (17.0 cm)

Internal volume: Approximately 350 millilitres.

Leakage: Less than 0.01 L/min at 3 psi (211 cmH_2O) internal pressure.

Construction: Filter housing and end cap material- Polysulfone plastic.

Connectors: ISO 22 mm male conical connectors.

WARNING

1. To ensure that all breathing circuit connections are leak-tight, complete a breathing circuit leak test according to ventilator operating instructions each time the Re/X700 filter is installed on a ventilator.

CAUTION

1. To ensure proper installation, observe the flow direction markings on the filter. Handle the filter with care; abusive or careless handling can damage the filter media. Discard any filter that rattles or has bulk liquid or debris visible inside the filter shell (visible bulk liquid may indicate breakdown of the hydrophobic barrier), or has visible rips or tears. Sterilize before disposal.

Note: Each Re/X700 filter is supplied nonsterile, and should be sterilized before use.

Cleaning and Sterilization

If it is necessary to clean the filter housing, wipe it down with a damp cloth. Any cleaning solution used for external wipe-down must be compatible with polysulfone plastic (filter housing and end cap material). The reusable Re/X700 filter is designed for sterilization by steam autoclaving at temperatures not to exceed 135°C (275°F). Effective sterilization occurs at a minimum cycle of 132°C (270°F) for 20 minutes for gravity-displacement cycles, or 132°C (270°F) for 4 minutes for pre-vacuum cycles (requires minimum 20 minute drying cycle after sterilization). Wrap the filter in muslin or equivalent paper before autoclaving. Follow the steam sterilizer manufacturer's instructions.

WARNING

1. **To avoid damage to the filter, do not immerse in liquid, or attempt to scrub or touch the filter media. Soaking or rinsing of the filter may cause retention of harmful residues and/or increased air flow resistance. Damage to the filter media can reduce its filtration efficiency.**
2. **Do not use ethylene oxide (ETO) to sterilize this filter. When handling the filter, always follow your hospital's infection control guidelines for handling infectious material.**

CAUTION

1. If nebulized medication is used, monitor filters more frequently as nebulized medication build up can occlude filter medium.


Note: Follow local governing ordinances and recycling plans regarding disposal or recycling.

Note: Covidien recognizes that cleaning, sterilization, sanitation, and disinfection practices vary widely among health care institutions. It is not possible for Covidien to specify or require specific practices that will meet all needs, or to be responsible for the effectiveness of cleaning, sterilization, and other practices carried out in the patient care setting. We do recommend that users of this product that require cleaning and sterilization/disinfection consider the Standards and Recommended Practices for sterilization published by Association for the Advancement of Medical Instruments (AMI) as well as the Center for Disease Control (CDC) publication: A guideline for Prevention of Nosocomial Pneumonia.

Filter Maintenance

Remove and steam sterilize the Re/X700 filter after each patient use. Look for any visible damage to the filter housing or media. Replace the filter if the filter body or the end caps are cracked, or if the filter media is visibly damaged. Check the filter flow resistance before each reuse. See Figure 1.

Figure 1: Pressure Drop Test Setup

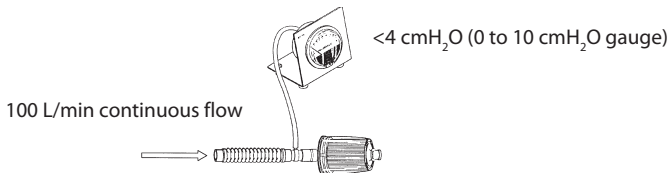Filter Replacement

The Re/X700 filter is designed to withstand up to 100 repeated autoclave sterilization cycles at temperatures not exceeding 135°C (275°F). Covidien recommends inspecting the filter before each reuse and discarding it if the filter housing or media is visibly damaged or if the filter flow resistance is too high for your application. Covidien recommends that use of the filter be discontinued if flow resistance is greater than 4 cmH₂O at 100 L/min flow.

Covidien recommends replacing the filter after a maximum of one year of service or after 100 steam autoclave cycles, whichever occurs first. For your convenience, a space is provided on the side of the filter to note anticipated retirement date. Calculate an expected filter retirement date based on your typical use rates.

Ordering Information

Item	REF	Description
1	G-060525-00	Re/X700 Reusable Expiratory Filter.
2	10076942B	Directions for use, Re/X700 Reusable Expiratory Filter

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© 2011 Covidien.
COVIDIEN, COVIDIEN with logo and Covidien logo are U.S. and/or internationally registered trademarks of Covidien AG. Other brands are trademarks of a Covidien company.
Covidien Ireland Limited,